

ELATION™

BALLOON DILATOR



EXPANDING THE POSSIBILITIES™



ELATION™ BALLOON DILATOR

Merit Medical Endotek's Elation balloon dilators provide precise, controlled, three stage dilation and crystal clear visualization during balloon endoscopy. Elation is designed for accurate sizing at three specific pressures, rapid and complete deflation as well as easy catheter withdrawal from the scope.

TAPERED, ATRAUMATIC TIP
Facilitates balloon placement in tight strictures while minimizing perforation risk

RADIOPAQUE MARKER BANDS
Allow for fluoroscopic visualization during balloon placement

PROPRIETARY CATHETER DESIGN
Provides improved pushability and rapid deflation

UNIQUE BALLOON DESIGN
Ensures three specific diameters at three corresponding pressures, providing accurate and repeatable balloon sizing

UNIFORM BALLOON EXPANSION
Consistent expansion over the full length of the balloon to reduce balloon waisting




Coupled with Endotek's BIG60™ inflation device, our wide array of Elation balloon dilators help alleviate common issues associated with balloon dilation.



Catheter securement clip allows for control of the catheter when not in use.



Marker bands indicate when the balloon dilator is ready to exit the scope and allow for an additional method for stricture measurement.

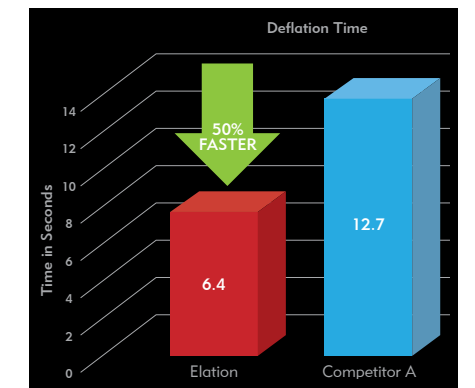

Improved deflation minimizes the risk of incomplete deflation and premature removal from the scope.*

*Data on file

RELATION SINGLE-USE FIXED WIRE BALLOON DILATOR

DISPOSABLE

Specifications (Packaged 1 unit)

Catalog Number*	Inflated O.D. (mm)	Balloon Length (cm)	Inflation Pressure (ATM)	Min. Working Channel (mm)	Catheter Size (F)	Working Length (cm)
E6	6-7-8	8	6 - 8.5 - 10	2.8	7.5	180
E8	8-9-10	8	4 - 6 - 8	2.8	7.5	180
E10	10-11-12	8	4.5 - 6.5 - 8	2.8	7.5	180
E12	12-13.5-15	8	4 - 6.5 - 8	2.8	7.5	180
E15	15-16.5-18	8	4 - 5.5 - 7	2.8	7.5	180
E18	18-19-20	8	3.5 - 5 - 6	2.8	7.5	180



Fixed Wire Balloon Dilator

* smart number
E=esophageal indication

RELATION SINGLE-USE WIREGUIDED BALLOON DILATOR

DISPOSABLE

Specifications (Packaged 1 unit)

Catalog Number*	Inflated O.D. (mm)	Balloon Length (cm)	Inflation Pressure (ATM)	Min. Working Channel (mm)	Catheter Size (F)	Working Length (cm)
EPB6	6-7-8	5.5	6 - 8.5 - 10	2.8	7.5	180
EPB8	8-9-10	5.5	4 - 6 - 8	2.8	7.5	180
EPB10	10-11-12	5.5	4.5 - 6.5 - 8	2.8	7.5	180
EPB12	12-13.5-15	5.5	4 - 6.5 - 8	2.8	7.5	180
EPB15	15-16.5-18	5.5	4 - 5.5 - 7	2.8	7.5	180
EB18	18-19-20	5.5	3.5 - 5 - 6	2.8	7.5	180
EPCB6	6-7-8	5.5	6 - 8.5 - 10	2.8	7.5	240
EPCB8	8-9-10	5.5	4 - 6 - 8	2.8	7.5	240
EPCB10	10-11-12	5.5	4.5 - 6.5 - 8	2.8	7.5	240
EPCB12	12-13.5-15	5.5	4 - 6.5 - 8	2.8	7.5	240
EPCB15	15-16.5-18	5.5	4 - 5.5 - 7	2.8	7.5	240
ECB18	18-19-20	5.5	3.5 - 5 - 6	2.8	7.5	240



Wireguided Balloon Dilator

* smart number
E= esophageal indication
P= pyloric indication
B= biliary indication
C= colonic indication

BIG 60 INFLATION DEVICE

Catalog Number	Description
ENDO-AN6012	60ml Inflation Device with Angled Gauge (Packaged 5 Per Box)



BIG60 INFLATION DEVICE
The all-in-one, easy-to-use design of our 60ml inflation device facilitates complete and rapid deflation, reducing procedure time and expense.

MAXXWIRE® GUIDE WIRE

Specifications

Catalog Number	O.D. (in)	Balloon Length	Shaft Tip Style	Packaged
80129-404	0.035	180	Super Stiff Straight	1 unit
80129-406	0.035	260	Super Stiff Straight	1 unit



www.merit.com/endotek

Merit Medical Systems, Inc. • 1600 West Merit Parkway • South Jordan, Utah 84095 • 1-801-253-1600 • 1-800-35-MERIT endotek.merit.com
Merit Medical EUROPE, MIDDLE EAST AND AFRICA (EMEA) • Amerikalaan 42, 6199 AE Maastricht-Airport • The Netherlands • Tel: +31 43 358 82 40
Free phone for specific country:
Austria 0800 295374 • Belgium 0800 729 06 (Dutch) 0800 731 72 (French) • Denmark 80 88 00 24 • France 0800 91 60 30 • Germany 0800 182 08 71
Republic of Ireland 1800 553 163 • Luxembourg 8002 25 22 • Netherlands 0800 022 8184 • Sweden 020 79 24 45 • UK 0800 973 115 • Finland 00800 6374 8633