


# ALIMAXX-B™

Biliary Stent System



EXPANDING THE POSSIBILITIES

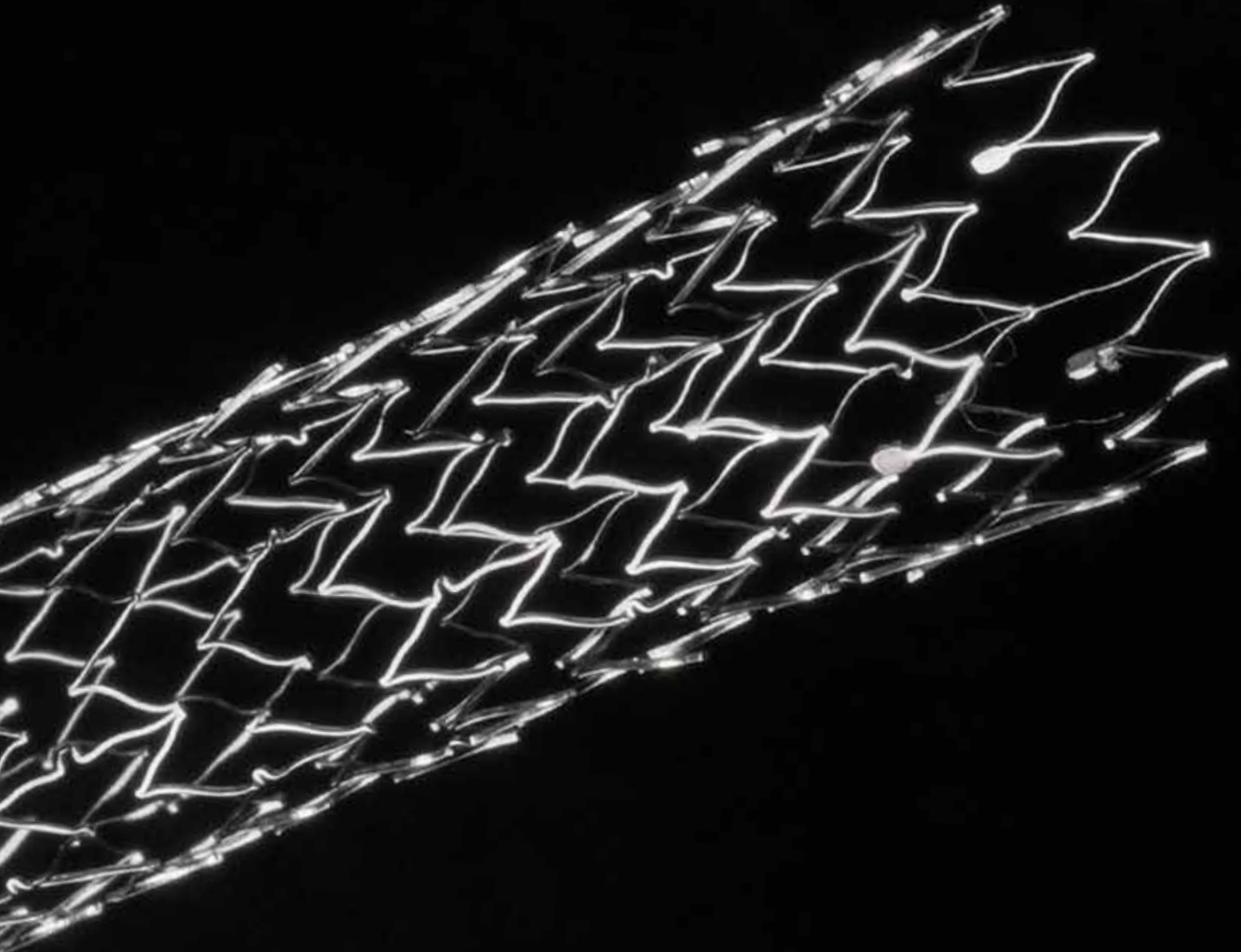


**MERTMEDICA**  
**ENDOTEK™**

# ALIMAXX-B™

## BILIARY STENT SYSTEM

*The uniquely designed metal stent, with its easy to use deployment system expands the possibilities in the treatment of malignant biliary strictures.*

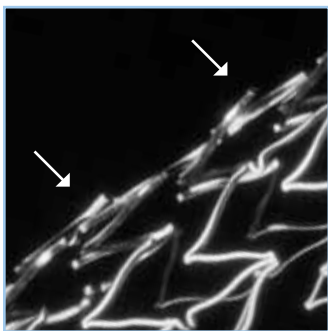


### **ANTI-MIGRATION STRUTS**


*Unique feature reduces stent migration.*



**LASER CUT NITINOL CONSTRUCTION** provides for virtually no foreshortening or elongation.

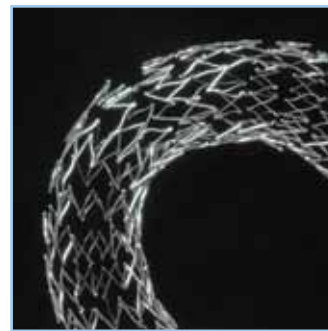

**ANATOMY SPECIFIC** High conformability aids in stent placement of tortuous anatomy and maintains luminal patency.

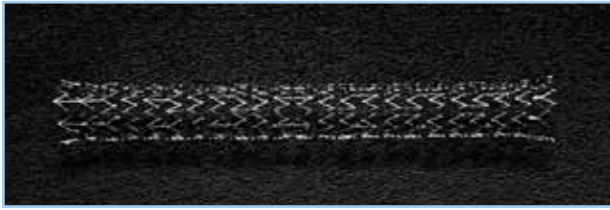


### **ANATOMY SPECIFIC GEOMETRIC SHAPE**

*The laser cut nitinol stent was designed to meet specific anatomical requirements of the biliary tree.*


*An open cell design allows for Y stenting at the hilar region.*

*The ALIMAXX-B's geometric shape provides increased and uniform radial force with true dynamic expansion.*

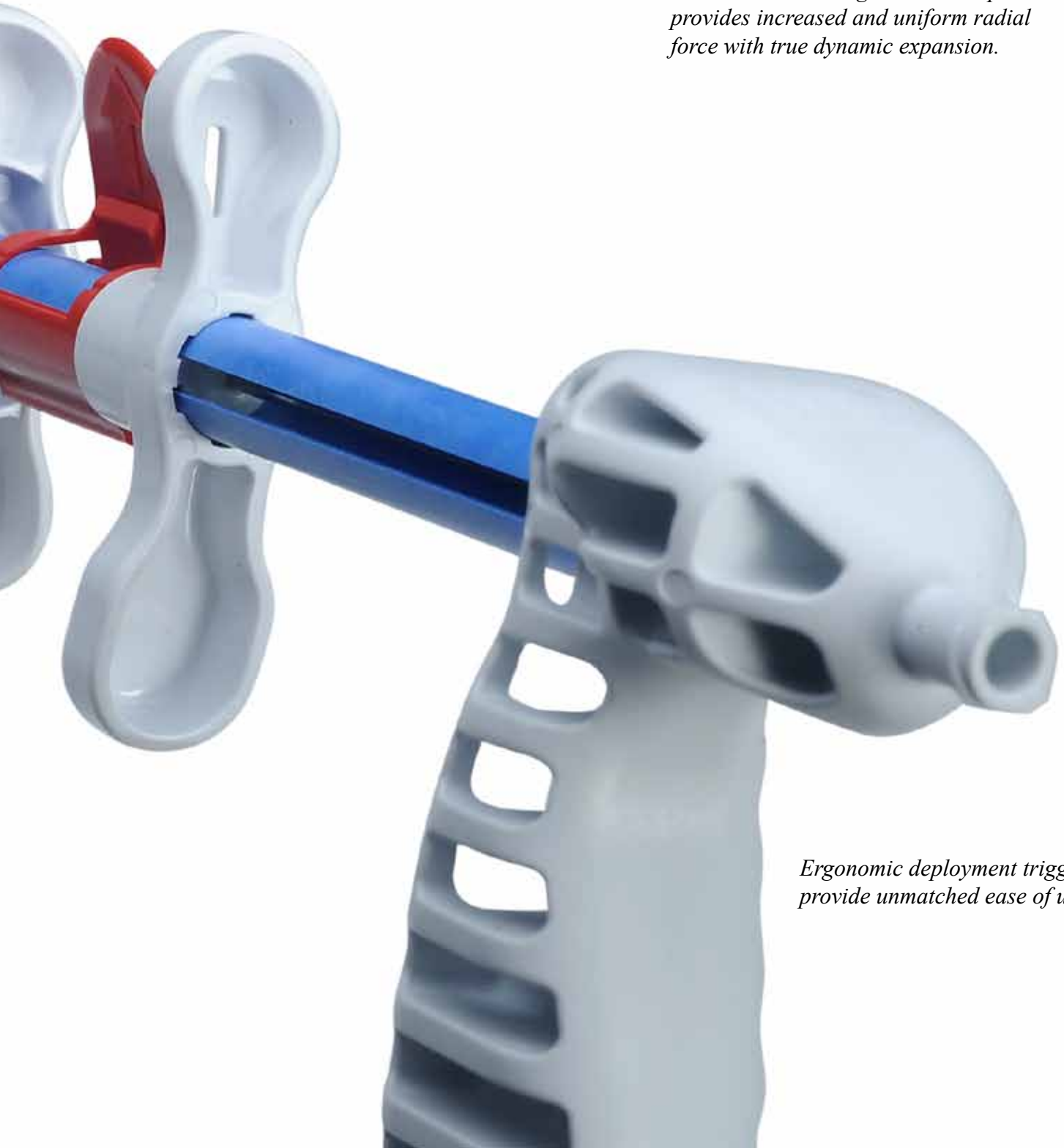


*Ergonomic deployment triggers provide unmatched ease of use.*

## FLUOROSCOPIC VIEW

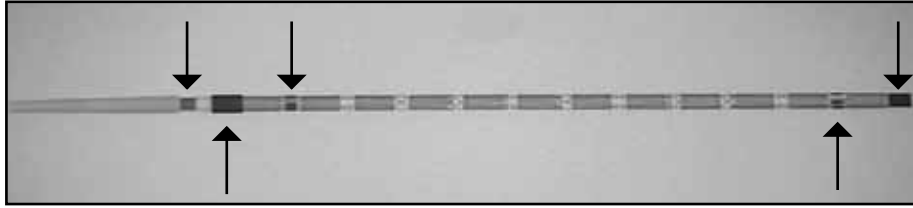

- 5 Radiopaque Markers
- One of two distal radiopaque markers moves with the outer sheath during deployment to help identify the progress of deployment.

## ENDOSCOPIC VIEW

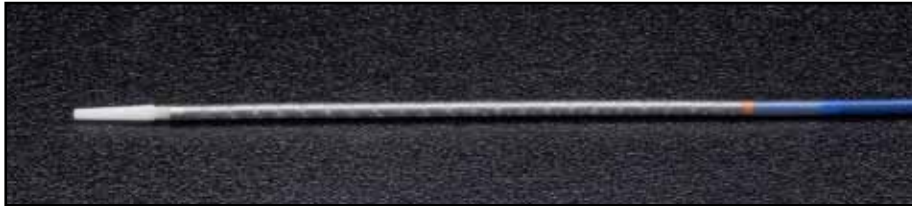

- Clear outer catheter allows increased visibility of proximal end of stent during deployment.
- Proximal end of stent is confirmed by orange endoscopic marker.

## ALIMAXX-B Metal Biliary Stent

Catalog Number	Stent Diameter	Stent Length	Delivery System Length	Delivery System Diameter
<b>70129-204</b>	8mm	40mm	190cm	6.5F
<b>70129-205</b>	8mm	60mm	192cm	6.5F
<b>70129-206</b>	8mm	80mm	194cm	6.5F
<b>70129-208</b>	10mm	40mm	190cm	6.5F
<b>70129-209</b>	10mm	60mm	192cm	6.5F
<b>70129-210</b>	10mm	80mm	194cm	6.5F


## ALSO AVAILABLE



Fully Covered  
ALIMAXX-ES™  
Esophageal Stent



MAXXWIRE™  
Guide Wire

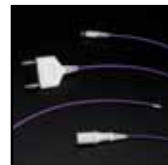

Brighton™ Bipolar  
Probe



ENDOTEK BIG60™  
Inflation Device



Merit Medical Systems, Inc. • 1600 West Merit Parkway • South Jordan, Utah 84095 • 1-801-253-1600 • 1-800-35-MERIT

[www.merit.com](http://www.merit.com)

Merit Medical EUROPE, MIDDLE EAST AND AFRICA (EMEA) • Amerikalaan 42, 6199 AE Maastricht-Airport • The Netherlands • Tel: +31 43 358 82 40

Free phone for specific country:

Austria 0800 295374 • Belgium 0800 729 06 (Dutch) 0800 731 72 (French) • Denmark 80 88 00 24 • France 0800 91 60 30 • Germany 0800 182 08 71

Ireland (republic) 1800 553 163 • Luxembourg 8002 25 22 • Netherlands 0800 022 8184 • Sweden 020 79 24 45 • UK 0800 973 115 • Finland 00800 6374 8633

401991001/A ID 111710